

The world's leading massage chair, Fujiiryoki.
Our technology has been evolving over time since 1954.

**Fujiiryoki, the world's first! mass-produced
massage chair, recognized as mechanical heritage!**

Certified August 7, 2014



Mechanical Heritage No. 69
Fuji Automatic Massage Machine,
The world's first mass-produced massage chair

FUJIIRYOKI
SINCE 1954 JAPAN

MASSAGE CHAIR

Made in Japan

JP-3000

CONDITIONING MASSAGE

**CYBER-RELAX JP-3000
Massage chair**

Specifications

- Dimensions:
(with tucked away footrest) about 86cm width,
141cm length, 124cm height
(reclined and with extended footrest) about 86cm width,
212cm length, 73cm height
(width including remote control holder: about 90cm)
- Minimum width before assembly: about 70cm
- Reclining angle: about 117~164 degrees
- Weight: about 96kg
- Power source: 100-240V (50/60Hz)
- Power consumption: 125W
- Continuous use time: 30 minutes
(automatic timer 7/18/30 minutes)



Black Gray (BG)

Smoke Beige (SB)


White (WH)

Emplacement • Installation


- Recommended moving space
70cm for the body + 6cm or more
for carrying

● Installation

Space needed for reclining has to
be accounted for and the chair
should be installed on a plane
surface. At least 30cm should be
left for footrest extension and
50cm should be left for backrest
reclining.

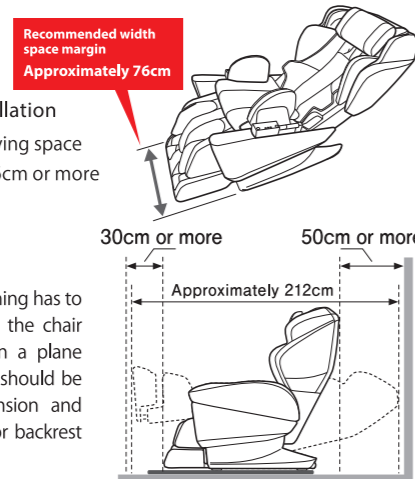

Dimensions (cm)

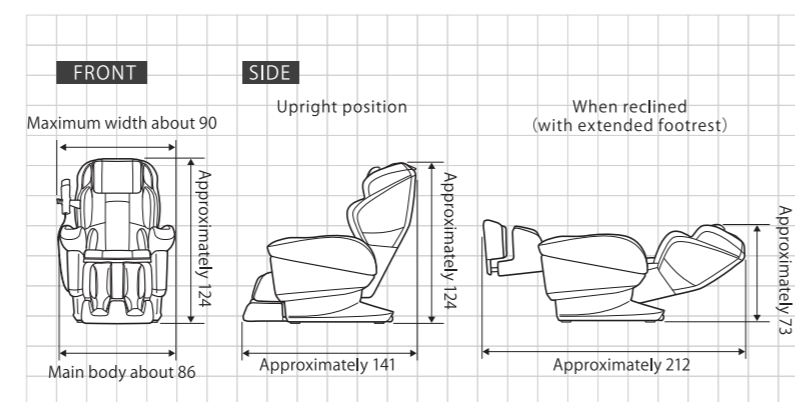

Purpose and Usage Effects

Hand massage substitute. For general home use.

- Allaying tiredness ● Promoting blood circulation ● Eliminating muscle fatigue ● Relieving muscle stiffness ● Relieving neuralgia and muscle soreness

- Do not use if you have the following conditions, otherwise it may cause abnormalities in your body.
- People who are not allowed by doctors to be massaged (examples: thrombosis (embolism), severe aneurysm, acute venous tumor, various dermatitis and skin infections (including inflammation of subcutaneous tissue), etc.)
- If you have any of the following conditions, please be sure to consult your doctor before use, otherwise it may cause abnormalities in your body.
- People who use devices susceptible to electromagnetic waves such as pacemakers or have implanted medical electronic devices in their bodies
- people with malignant tumors ● people with heart disease ● women who are pregnant or have just given birth ● people with high peripheral circulation impairment due to diabetes ● people with skin trauma ● people who need to keep their bodies and minds quiet ● people with a body temperature of 38°C or higher (fever) (e.g., people with acute inflammatory diseases [malaise, chills, blood pressure changes, etc.] and who are in a period of intense symptoms), people with osteoporosis, spinal fractures, acute [painful] diseases, ● people with abnormalities of the backbone (spine) or scoliosis ● people with sprains, strains and other inflammatory symptoms ● people with herniated discs ● people with tenosynovitis ● people who are considered to have a loss of temperature sensation ● people with other special physical symptoms, or people who are being treated at medical institutions
- Do not allow it to be used by children alone or by people unable to express their own will. Please also keep children away.

The Health and Beauty Manufacturer
FUJI MEDICAL INSTRUMENTS MFG.CO.,LTD.

14F, Oe Building, 1-22, Noninbashi 1-chome,
Chuo-ku, Osaka 540-0011, Japan
www.fujiiryoki.com/en



サイバーリラックス

CYBER-RELAX

FUJI MEDICAL INSTRUMENTS MFG.CO.,LTD.

●The design and specifications are subject to change without notice due to product improvement and other factors. ●Please note that there may be color differences between the print and the actual product. ●Photos and illustrations are for reference only. ●The contents of this catalogue are as of March 2023.

A whole new massage experience that allows both body and mind to be at their best.

There have been dramatic changes in the way we live and work all over the world. Various types of fatigue and environmental changes can also cause mental stress. In order to be active every day and be at one's very best, it is necessary to take care of both body and mind.

CONDITIONING MASSAGE

A conditioning massage created to take care and tune up physical, mental and health state.




1
Massage care through AI analysis

BODYCARE

Using AI analysis to run the massage that best suits the body's current state.

- AI Double Sensing +
- 5D-AI Mecha PLUS
- Foot sole Shiatsu
- 38 Course Menus
- User Customize Mode

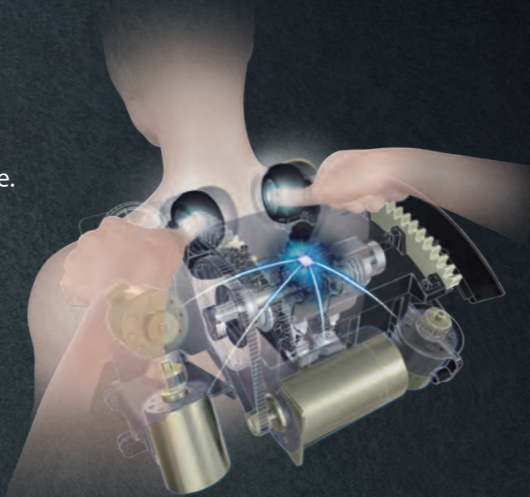


2
Leading to the best state

BALANCE

Using a variety of stretches to adjust and balance the muscles.

- Air Massage
- 5 Full Body Stretch Menus
- 6 Body Parts focused Stretch Menus

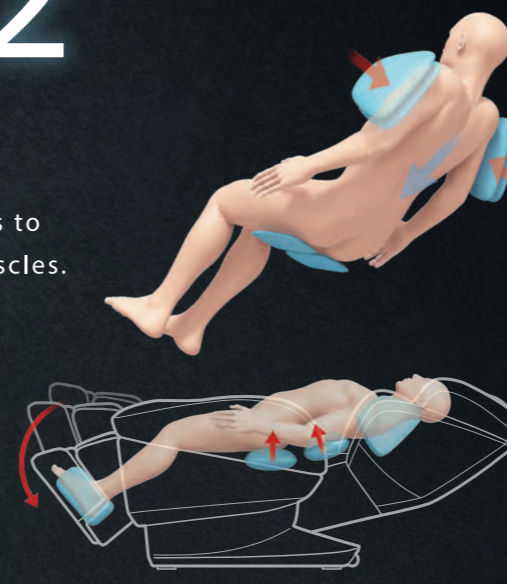

3
A relaxing experience for both Body and Mind

BEING


Creating moments of relax to take care of the mental state.

- Mindfulness Course
- Total Relax Position
- Double Heating
- Bluetooth® Supporting Speakers



BODYCARE

To obtain the best bodycare, the accurate AI sensing technology allows to optimize the massage according to each user.

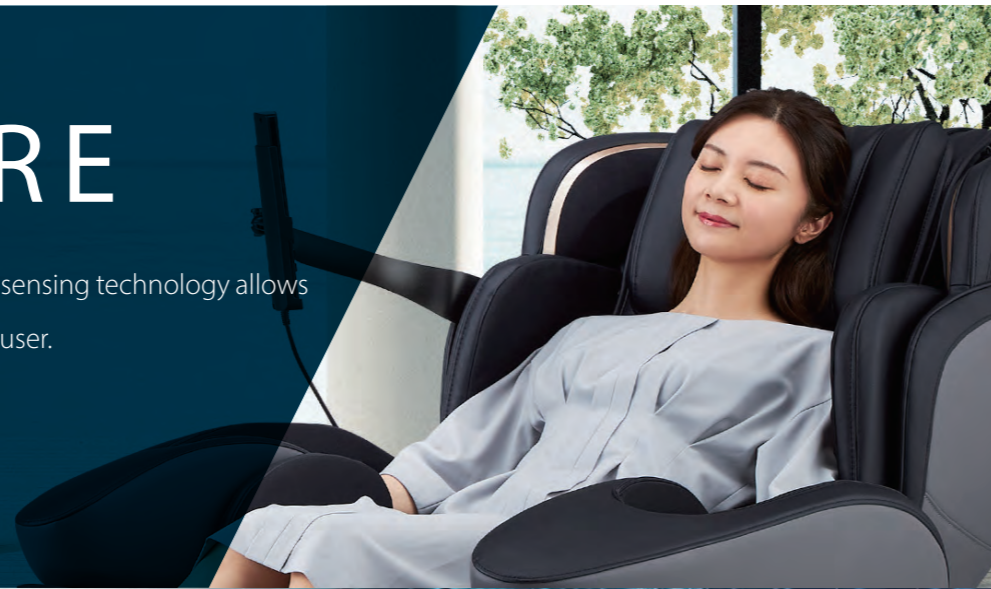

5D-AI PLUS




Fujiiryoki AI Technology

Through the AI, the control system automatically passes the user's body shape, muscles status and other information into the massage. Intensity, techniques, smoothness are adjusted and implemented according to each user.

AI By detecting 8 backline massage points, it's possible to automatically optimize and adjust the massage according to the user.

NEW AI Double Sensing +



The "Shoulder Position Detection System", which detects the position of the shoulders, calculates the massage point, then the "Backline Detection System" detects the backline by measuring 8 different points along the spine. The 5D-AI movement PLUS automatically adjusts itself along this spinal curve to achieve the most comfortable massage.

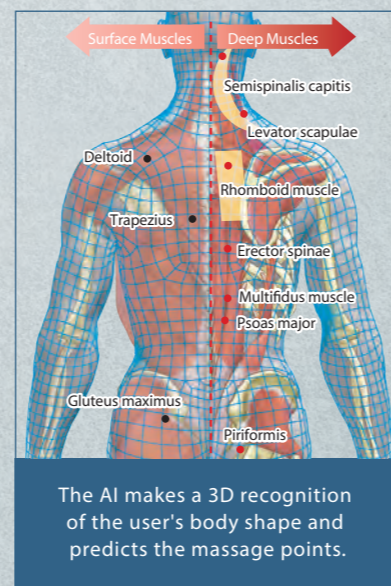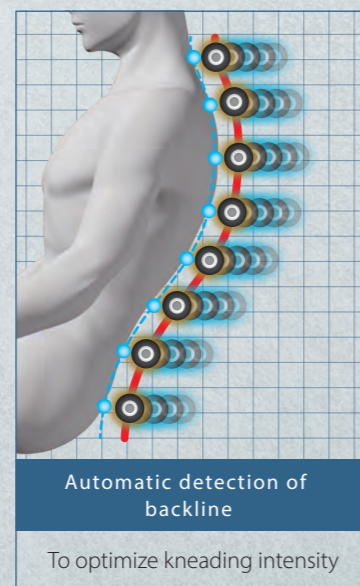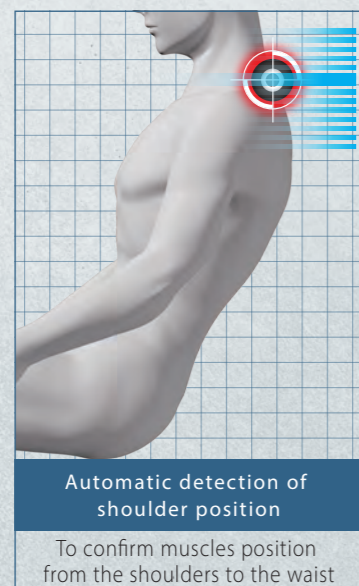


AI The two kneading balls mecha has been cultivated by Fujiiryoki for more than 20 years. Fine and efficient, Japanese quality kneading is achieved through the AI.

NEW 5D-AI Mecha PLUS

Through the AI control of kneading and a high-precision motor, the past mechanical movements have been successfully upgraded into a real hand-like, smooth massage. Moreover, intensity is adjusted according to the muscle state of each person in real time to relieve body stiffness with a stable and on-point massage.

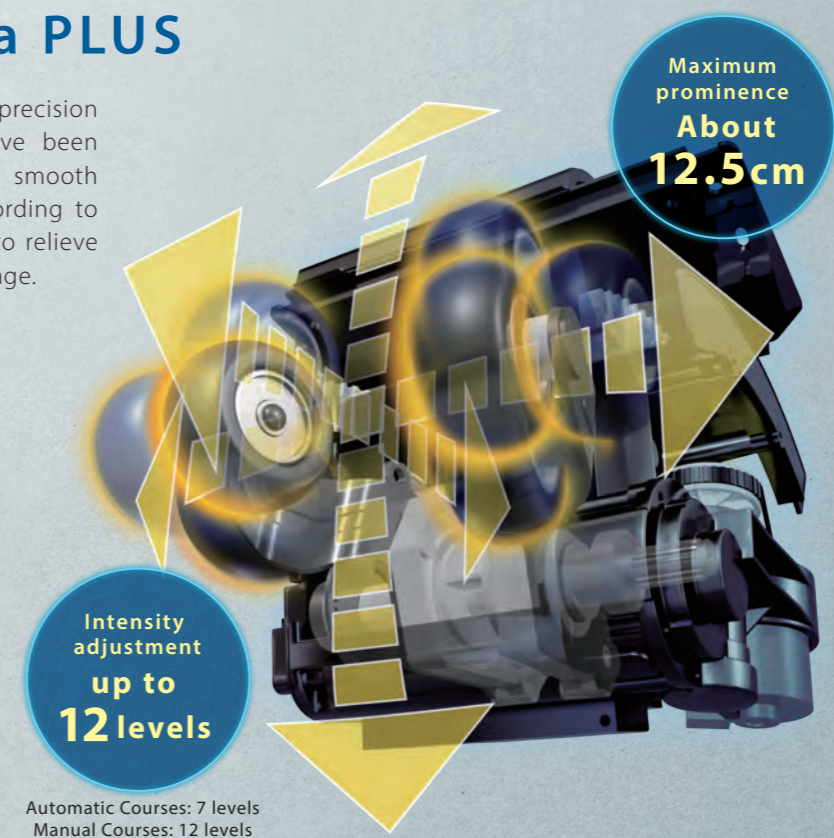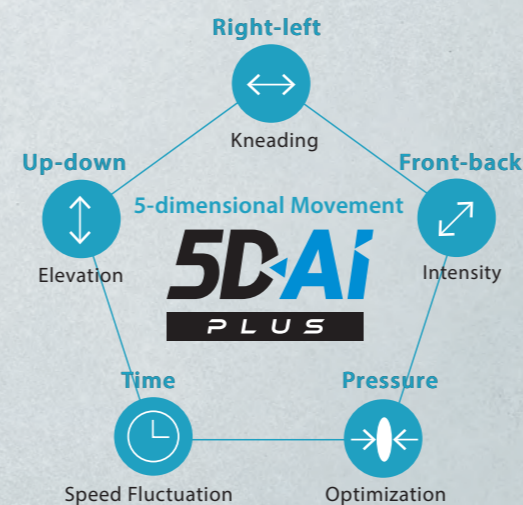


Kneading action of 5D-AI Mecha

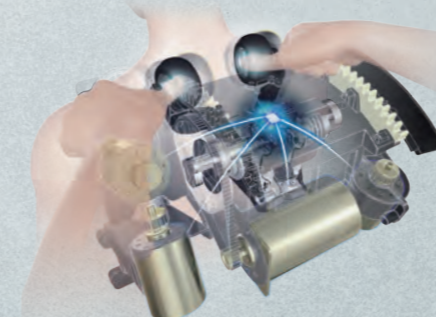

Smooth Approach



Thorough Shiatsu

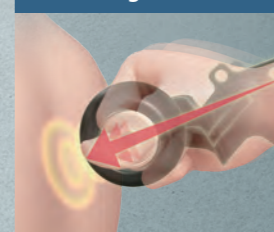


Firm Kneading




Foot Sole Shiatsu

The airbag keeps ankle and instep in place, as if a hand grips the foot, and then a special roller massages the sole of the foot.

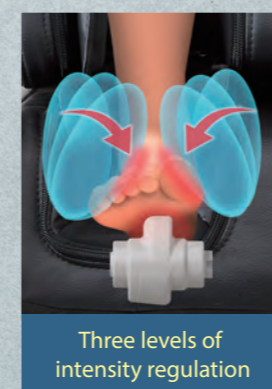

Foot Roller Massage

By regulating the strength of the air bag holding ankle and instep, the intensity of the foot roller can be adjusted in three levels. According to preferences, the foot roller can also be turned off.

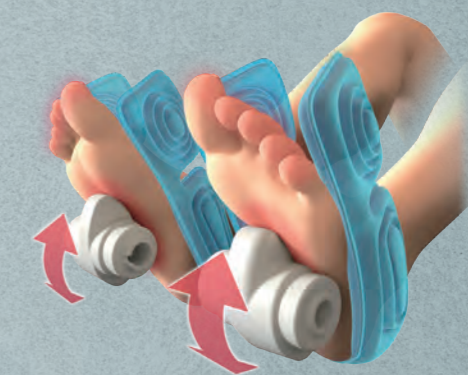

※"Shiatsu" is here used to describe the roller movements

38 Course Menus creating a special experience every day

A rich offer to satisfy all preferences, as it's possible to choose among 38 courses, 23 automatic and 15 part-focused ones.

The large screen is easy to read and operate

4.3" Full Color Touch Screen Remote Control

An easy-to-understand voice guidance helps with each operation step.

Voice Guidance Function

The impressive number of menus allows to choose the massage course according to daily fatigue levels.

23 Types of Automatic Courses

30-minute Courses

| | | | | |
|---|--|--|---|---|
| VIP Course Concentration on body parts prone to stiffness. Providing a salon-like time of premium relaxation. | Elite Therapy Enjoy a luxurious and comfortable massage as if you were being massaged by multiple masseurs at the same time. | Night Healing An effective stretch and massage course, ideal before bedtime. | Conditioning NEW A stretch course focused on slowly and thoroughly relaxing muscles of the whole body. | Mindfulness NEW A slow and comfortable Air massage to eliminate fatigue and tiredness of the whole body. |
|---|--|--|---|---|

18-minute Courses

| | | | | |
|---|---|---|--|---|
| Whole Body Course Extreme Rejuvenation Rejuvenation Refresh Soft | Neck/Shoulder Course Extreme Rejuvenation Rejuvenation Refresh Relax | Waist Course Extreme Rejuvenation Rejuvenation Refresh Relax | Stretch Course Extreme Whole Body Stretch Whole Body Stretch Whole Body Air Technical Stretch | 7-minute Course Quick Course Course effectively relieving the whole body in a short time. |
|---|---|---|--|---|

7-minute Short Courses that focus on soothing specific body parts

15 Part Focus Techniques Unique Functions

6 Kneading Massage Techniques

| | | | |
|------------------------------|-----------------------|---------------|------------------|
| Neck Focus Massage | Waist Focus Massage | Neck Stretch | Shoulder Stretch |
| Shoulder Focus Massage | Sciatic Focus Massage | Back Stretch | Leg/Knee Stretch |
| Shoulder-blade Focus Massage | Foot Focus Massage | Waist Stretch | Arch Stretch |

3 Body Trunk Techniques

| | | |
|-----------------------|------------------------|---------------------|
| Back/Body Trunk Relax | Waist/Body Trunk Relax | Pelvis Core Stretch |
|-----------------------|------------------------|---------------------|

The "power" and "technique" of the two kneading balls mecha allow a precise massage.

53 Kneading Massage Techniques



Wavelet
Kneading and tapping at the same time with fine movements, massaging a wide area.



Tapping Spiral Knead **NEW**
Using small circular motions to knead and soothe stiff muscles, then tapping to loosen them.

Rubbing Massage **NEW**
Kneading while rubbing, smoothly relaxing the muscles.

Irregular Waves Wavelet **NEW**
Deep muscle massage through quick kneading and tapping.

Tapping Precision Spiral Knead **NEW**
Small circular motions to relax stiff muscles.

Upward Lifting/Downward Pressing **NEW**
Kneading slowly and deeply while stretching the muscles through pressing and lifting movements.

A customized massage chair to suit each user preferences

NEW User Customized Mode


User's shoulders position, back line, airbags intensity and reclining angle can all be registered.



step1 Body Shape Information Detection

AI Double Sensing +

Thanks to the "Shoulder Position Detection System", which detects the position of the shoulders, and the "S-line Detection System", which detects the back line, the AI is able to optimize the massage points automatically.

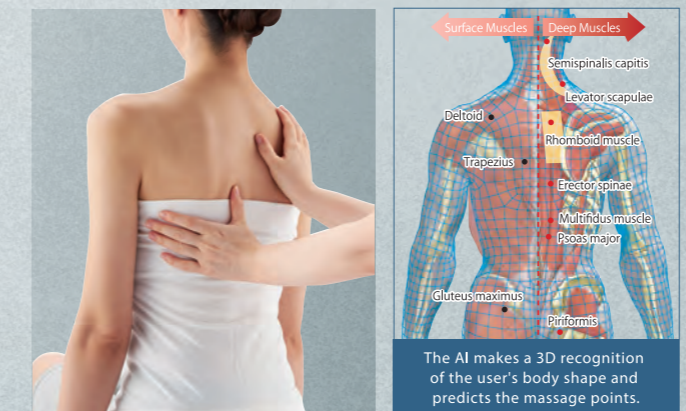

step2 Registration of Kneading Intensity from Shoulders to Hips

Customized Back Line

Based on the detected back line, a uniform massage is applied from the shoulders to the buttocks. In addition, the kneading intensity of each of the 8 back line points can be adjusted according to preferences.



step3 Registration of Most Relaxing Position

Customized Reclining Angle


Any position deemed to be comfortable by the user will be reflected and registered into all programs. The reclining angle can be adjusted according to preferences.



step4 Adjustment of Air Massage Intensity

Customized Air Massage

The intensity of each of the airbags that massage 8 different parts of the body can be registered independently, in real time.

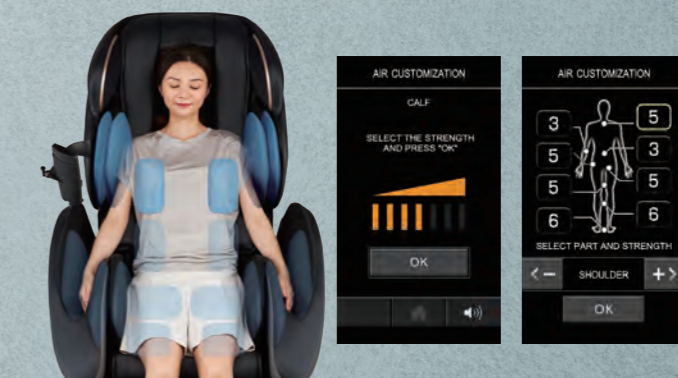

BALANCE


A massage that relaxes and stretches the muscles, creating a body that is balanced and less likely to feel fatigued.




Through a carefully calculated airbag configuration and a fine strength adjustment function, the air massage perfectly fits the whole body.

High-efficiency Air Pressure System

The strength of the airbags can be adjusted in 7 levels for 8 parts of the body. It's possible to experience the Pulse Mode, that with very regular action promotes blood circulation, as well as the Air Massage Mode, that through different speeds imitates a real hand massage.

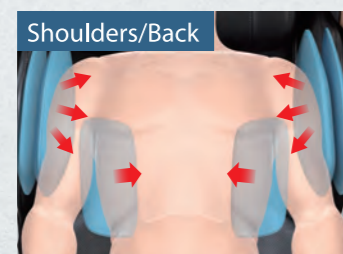

Shoulders/back Massage
(7 levels of adjustment)
A wrapping sensation massage, from the shoulders to the back.

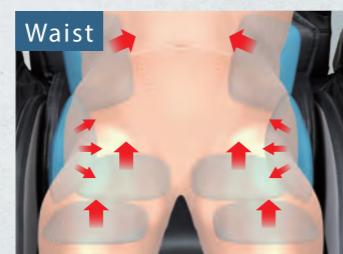

Waist Massage
(7 levels of adjustment)
Loosens a wide range of muscles, from around the pelvis to the back of the knees.



Arm Massage
(7 levels of adjustment)
Massage from the elbow to the fingertips.



Legs Massage
(7 levels of adjustment)
Through repetitive compression and release, gives relief to calves and ankles.



Exclusive Function

High-speed Inflation Airbag System

The built-in high-flow and high-pressure air pump allows the high-capacity airbags to operate powerfully, enhancing the massage effect.

Stretching and relaxing muscle tension of the whole body and of single body parts for a perfect tune-up.

Diverse Stretching Menus

5 Types of Full Body Stretch

Stretching that carefully tunes up the whole body (30 minutes)

- Conditioning Course

Everyday full-body Stretch Menu (18 minutes)

- Technical Stretch
- Extreme Whole Body Stretch
- Whole Body Stretch
- Whole Body Air




6 Types of Region-focused Stretch

1 Neck Stretch

2 Shoulder Stretch **NEW**

Stretch menu for shoulders, obtained through lifting up and releasing movements.

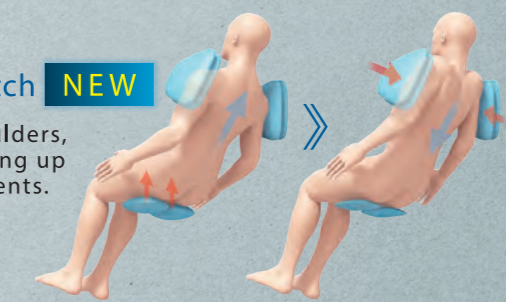


3 Back Stretch

4 Waist Stretch


5 Arch Stretch **NEW**

Menu that slowly bends and stretches the upper body

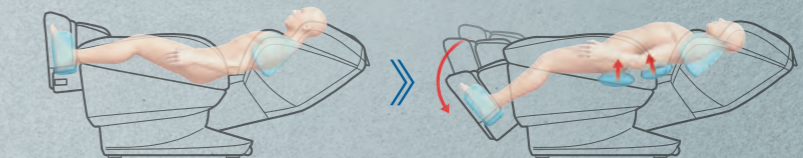

6 Leg/Knee Stretch **NEW**

Stretch menu that focuses on the knees

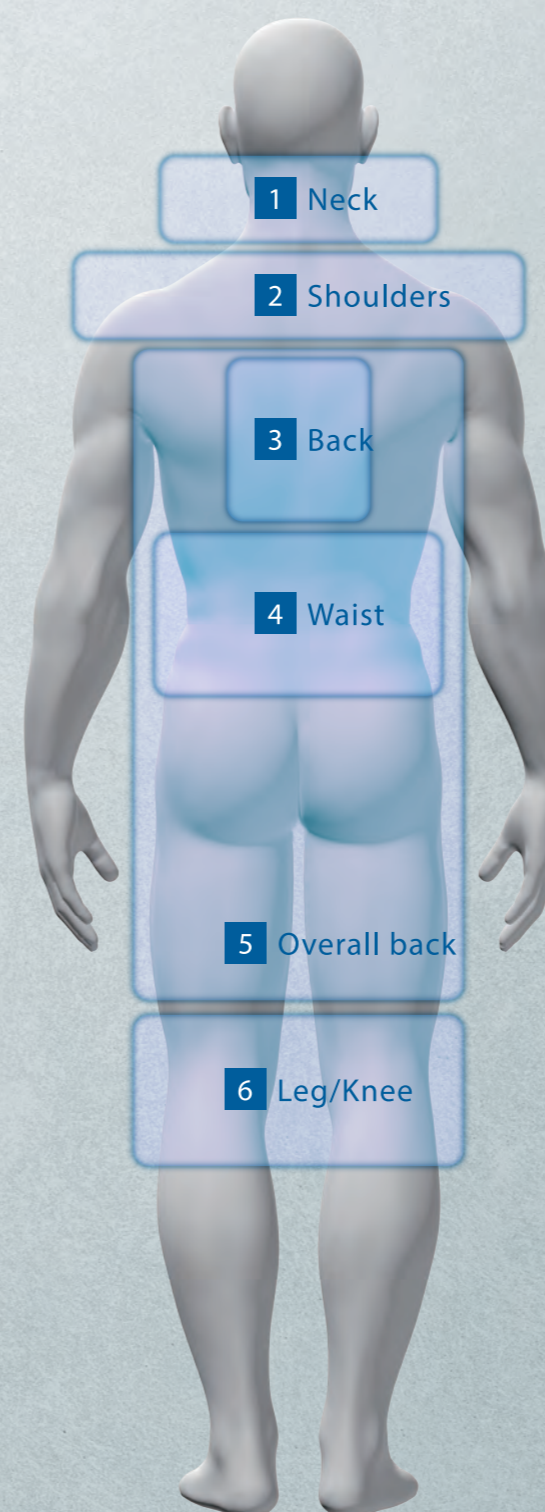

Power Stretch

Exclusive Function

Double Stretch

While pelvis and legs are kept in place through the airbags, the legs are also pulled forward and downward, stretching the body and relieving muscles stiffness.



Exclusive Functions

Knee Stretch

Legs Stretch

BEING

Creating moments of relaxation and enjoyment to bring tranquility to the mind and guide the body to its best-performing state.



Bluetooth® Supporting Speakers

Specially designed left and right bass-reflex speakers are placed near the ears to achieve a luxurious sound-surrounded space.

Clear Voice Sounds

Built-in Multi-band Equalizer

Stereo Sound

Built-in DSP Digital Speakers

Rich Output Sound Range

Full-range Speakers


The Bluetooth® typeface and logo are registered trademarks of Bluetooth SIG, Inc. and Fujiiryoki Ltd. is authorized to use the above marks. All other trademarks and company names are the property of their respective owners.

Double Heating

Two heating systems are built in to create a more comfortable massage time.

Foot Sole Heating

Gently warming up the feet, which tend to get cold easily.

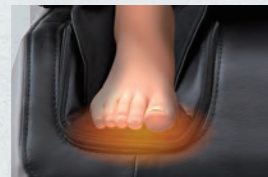

*The warming function can be used independently.
*Set temperature: about 43 degrees (set temperature is different from perceived temperature.)

Exclusive Function

Quick-Heat Multi-heaters

The pads with built-in warming function can be placed behind the back to warm up around the waist, or can be unzipped and placed on the front to warm around the abdomen.



Waist warming



Abdomen warming



NEW Mindfulness Course [automatic]

Based on the concept of "meditation", a 30-minute massage performed through airbags only. The slow air massage carefully and gently soothes the entire body.



NEW Total Relax Position

The relaxation level is enhanced by the comfortable lying position, also perfect to optimize the stretching effects. You can actually feel the stretching effect.



BLACK GRAY
(BG)



SMOKE BEIGE
(SB)



WHITE
(WH)



The history of massage chairs is the history of Fujiiryoki.

The world's leading massage chair, Fujiiryoki. Our technology has been evolving over time since

1954



The "starting point" of massage chairs. Machine No. 1 "Fuji Automatic Massage Machine"

2001



CYBER-RELAX AS-001
Built-in 3D massage point detection system

2003



CYBER-RELAX AS-800
Built-in shoulder and foot sole massage functions

2015



CYBER-RELAX JP-1000
Built-in deep tissue massage mecha PRO (Kiwami)

2017



CYBER-RELAX JP-1100
Built-in foot sole Shiatsu

2019



CYBER-RELAX JP-2000
Built-in 5D-AI mecha

2023



CYBER-RELAX JP-3000
Built-in User Customized Mode